ESPAÑOL

Balanza Body Fat

PK-DS46

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PUNKTAL

WWW.PUNKTAL.COM.UY

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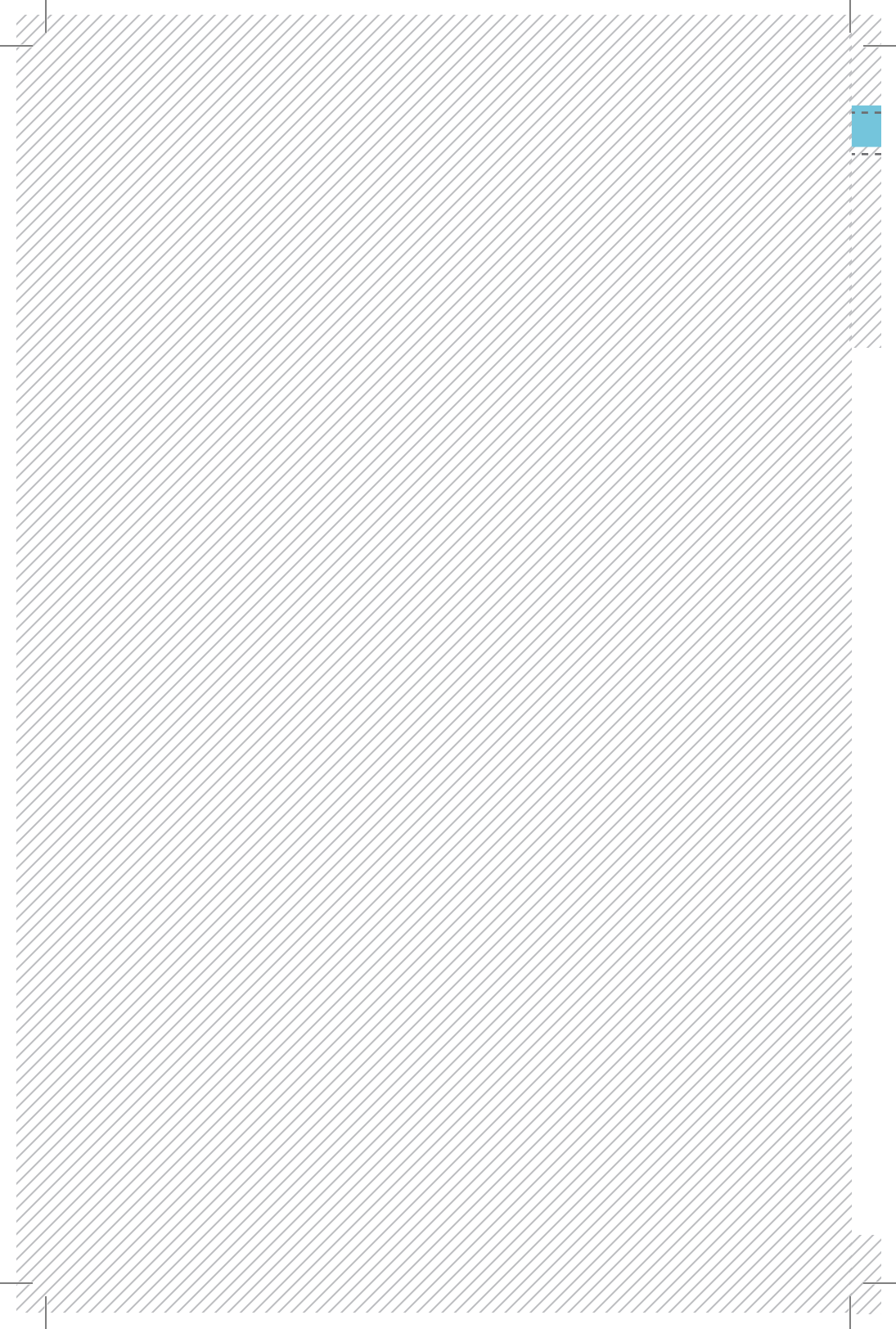
PUNKTAL

WWW.PUNKTAL.COM.UY

Thank you for choosing one of our products . We congratulate you for your purchase.

Our appliances are designed to meet the highest standards of quality.

GENERAL SAFETY INSTRUCTIONS:

Read the operating instructions carefully before putting the appliance into operation and keep the instructions including the warranty, the receipt and, if possible the box with the internal packing.

- The device is designed exclusively for private and not for commercial use.
- Only use the device for the intended purpose.
- Do not use a damaged device.
- Keep away from heat, direct sunlight, damp (never immerse in liquids) and sharp edges.
- Do not try to repair the appliance on your own. Always contact an authorized technician.
- Pay careful attention to the following "Special Safety Instructions".

I) SYMBOLS IN THESE INSTRUCTIONS FOR USE

Important information for your safety is specially marked. It is essential to comply with these instructions in order to avoid accidents and prevent damage to the machine:

WARNING

This warns you of dangers to your health and indicates possible injury risks.

⚠ CAUTION:

This refers to possible hazards to the machine or other objects.

ℹ NOTE:

This highlights tips and information

II) CHILDREN

In order to ensure your children's safety, please keep all packaging (plastic bags, boxes, polystyrene etc.) out of their reach.

CAUTION

Do not allow small children to play with the foil as there is a danger of suffocation!

SPECIAL SAFETY INSTRUCTIONS FOR THIS MACHINE

WARNING

If plastic surface is moist or covered with water droplets, please wipe the surface dry as there is a danger of slipping.

CAUTION

- Handle the balance with care, as it is a precise measuring instrument. Do not prop it or jump on it.
- Do not load the appliance with more than 180kg to avoid damaging the device.
- If the scales are carried or not going to be used, please store them only in a horizontal position.
- The body fat analysis scale is not suitable for the following persons:
 - Children under 10 years and adults over 85 years
 - Persons with symptoms of fever, oedema, or osteoporosis
 - Persons undergoing treatment for dialysis
 - Persons who use a pacemaker
 - Persons who take cardiovascular medication
 - Pregnant women
 - Sportspersons who perform over 5 hours of intensive sporting activities each day
 - Persons with a heartbeat lower than 60 beats/minute

I) HANDLING BATTERIES

WARNING

- Do not expose the battery to high temperatures or direct sunlight. Never throw batteries into the fire. There is a danger of explosion!
- Keep batteries away from children. They are not a toy!
- Do not open the batteries by force.
- Avoid contact with metallic objects (rings, nails, screws etcetera.) there is a danger of short-circuiting!
- As a result of a short-circuit batteries may heat up considerably or even catch fire. This may result in burns.
- For your safety the battery poles should be covered with adhesive strips during transport.
- If a battery leaks do not rub the liquid into contact with the liquid they should be washed thoroughly and eyes rinsed out with clear water, if symptoms persist please consult a doctor.

II) INSERTING THE BATTERIES

i NOTE

Where required, the supplied batteries have been protected for transport with plastic foil.

1. Remove this foil before inserting the batteries into the appliance
2. Open the battery compartment on the back.
3. Insert 4 batteries of type AAA/R03. Please ensure that the polarity is correct(see the loor of the battery compartment)!
4. Close the battery compartment again.
5. If the device is not used for a prolonged period of time,please take the batteries out in order to prevent leakage of battery acid.
6. If the message "LO"appears on the display,replace the batteries.

△ CAUTION

- Different battery types or new and used batteries must not be used together.
- Batteries are not to be disposed of together with domestic wastes. Dispose of the used batteries at the relevant collection bin or at your dealer's.

III) LOCATION OF CONTROLS

1. Contact surfaces
2. Display
3. DOWN button ▼
4. UP button ▲
5. SET button

Display Symbol	Meaning
(Kg/lb/st)	Body weight
P0-9	User memory space
	Female
	Male
CM (FT/In)	Height
AGE	Age
FAT	Body fat
TBW	Body water
MUS	Muscle mass
BONE	Bone weight
%	Percentage-body fat/body water/ muscle mass/bone weight

TECHNICAL SPECIFICATIONS

Model:	PK-DS46
Power supply:	Batteries 1.5V type AAA/R03
Net weight of device:	2 Kg
Memory capacity:	10
Mass alignment:	from 75 to 225 cm (2'5" - 7'4 1/2" - ft/in)
Age:	from 10 to 85 años
Minimum weight:	from 2.5 Kg
Maximum weight:	150 Kg
Weight measurement unit:	100 gr increments
Body fat measuring unit:	0.1 %
Body fat measuring range:	4.0 % - 50.0 %
Body water measuring range:	35.0 % - 85.0 %
Muscle mass measuring:	1.0 % - 50.0 %

Subject to technical changes without prior notice.

This device has been tested according to all relevant current CE guideliness, such as electromagnetic compatibility and low voltage directives, and has been constructed in accordance with the latest regulations.

INFORMATION ON MEASURING BODY FAT

CAUTION

Do not use the scales if you use a pacemaker or another medical implant with electronic components. The function of the implant may be impaired the measurement.

i NOTE

- The device can only function correctly if you stand on the contact surface with bare, dry feet.
- The main principle of this device is based on the measurement of the electric impedance within the human body. In order to determine the actual body date, the device sends a weak, unnoticeable electric signal through your body.
- This signal is generally safe! However, please read our special safety instructions!
- The method of measurement is referred to as "bioelectric impedance analysis" (BIA), a factor which is based on the relationship between body weight and associated body water percentages as well as on other biological data (age, gender, body weight).
- There are heavy persons with a well-balanced percentage between fat and muscle mass, and there are seemingly slim persons with significantly high body fat values. Fat is a factor that puts serious strain on the body, in particular the cardio-vascular system. Therefore, controlling both values (fat and weight) is a prerequisite for health and fitness.

BODY FAT:

Body fat is an important component of the body. From protecting the joints as a vitamin repository to regulating the body temperature, fat plays a vital role in the body. Thus the goal is not radically reduce the proportion of body fat, but rather to achieve and maintain a healthy balance between muscle and fat.

GUIDELINE VALUES FOR BODY FAT AND BODY WATER CONTENT AS A PERCENTAGE(APPROXIMATE FIGURES)

Age (Years)	Women		Men		
	Body fat	Body water	Body fat	Body water	
< = 30	4.0 – 20.5	66.0 – 54.7	4.0 – 15.5	66.0 – 58.1	Underweight
	20.6 – 25.0	54.6 – 51.6	15.6 – 20.0	58.0 – 55.0	Normal
	25.1 – 30.5	51.5 – 47.8	20.1 – 24.5	54.9 – 51.9	Overweight
	30.6 – 50.0	47.7 – 34.4	24.6 – 50.0	51.8 – 34.4	Obese
> 30	4.0 – 25.0	66.0 – 51.6	4.0 – 19.5	66.0 – 55.4	Underweight
	25.1 – 30.0	51.5 – 48.1	19.6 – 24.0	55.3 – 52.3	Normal
	30.1 – 35.0	48.0 – 44.7	24.1 – 28.5	52.2 – 49.2	Overweight
	35.1 – 50.0	44.6 – 34.4	28.6 – 50.0	49.1 – 34.4	Obese

Underweight mark: -

Normal Mark: o

Overweight Mark: +

Obese mark: ++

USE

The device has a very sensitive sensor.If you want to adjust the settings,simply touch the buttons gently with your finger.if the device does not respond to touch,clean the glass surface with a soft,dry cloth.

i NOTE

The scales are factory set with the units "kg" and "cm".

To change the settings to "lb" or "kg" and "Ft:in", proceed as follows:

On the back side of the scale beside the battery compartment a switch can be found, push the below switch, choose between the measuring units kg/lb/st.

WEIGHING FUNCTION

- Place the scale on a hard, flat surface(avoid carpets). An uneven floor contributes to inaccurate measurements.
- Stand on the scales with both feet. The device switches on automatically.
- Distribute your weight equally and wait for the scale to determine your weight.
- The digits on the display will flash before your exact weight will be shown.
- In order to extend the battery life, the device will automatically turn off 10 seconds after you have stepped off the scale.

SAVING PERSONAL DATA

i NOTE

- The scale has default values that you can change, as described below.
- Carry out the following steps in sequence and within a few seconds. Otherwise, the device will automatically cancel input mode.
- Press the button SET. The memory space will flash on the display (USER0). Confirm the setting by pressing SET or select another memory space by pressing the UP/DOWN button. Then press SET to go to the next setting.
- The symbol for "male" or "female" will flash. Confirm the setting by pressing SET or use the UP/DOWN button to choose the symbol for your gender. Then press to go the next setting.
- The height setting will flash. Select your height by pressing UP or DOWN button. Holding down one of the buttons will accelerate the process. Confirm your entry by pressing SET.
- The setting AGE will flash. Select your age by pressing UP or DOWN button. Holding down one of the buttons will accelerate the process. confirm your entry by pressing SET.
- The device can only function correctly if you stand on the device with bare, dry feet.
- If necessary wipe the contact surface improve results.
- Due to the significant fluctuation of body fat, it is different to attain a standard value. On average the optimal testing time is between 6 p.m. And 8 p.m, since the body fat remains very consistent during this period of time.
- Press the SET button. The preset stored position date appears on the display (USER flash). Use the UP/DOWN button to select the stored position at which your date is saved.

- Wait for the USER display and 0.0kg (0.0lb/st) to stop flashing.
- Stand on the contact surfaces(1) of the scales with both feet.
- Firstly your weight will flash on the display.
- After a few seconds, the following results will be shown:
 - the body fat content in%(FAT symbol)
 - the body water content in %(TBW symbol)
 - the results analysis in the form of a bar
 - the muscle mass content in %(MUS symbol)
 - the bone mass content in %(BON symbol)
 - the calorie requirement per day
 - the present gender
 - the present type
 - the present height
 - the present age(AGE)

DISPLAY ERROR MESSAGE

If an error message appears on the display, wait until the device is turned off. Turn it on again to restart the operation. If the error display has not disappeared, take the batteries out of the battery compartment for 10 seconds.

- Clean the device with a slightly moist cloth. make sure that no humidity or waters seeps into the device.
- Use a commercially sold detergent or soap for cleaning.
- Dry the device well after cleaning with a dry and soft cloth.

8888	_____	Self test (when switching on) to check that the numeric display functions properly las funciones del display.
----	_____	The scale stands faultily.
Lo	_____	Low battery.
Err	_____	The percentage of body fat content lies above or below the limits of the scale.
EEEE	_____	The weight lies above the allowed maximum weight.

PK - DS46

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Salida: RMS 3WTHD=10%

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FM: 87.S•108MHZ

PUNKTAL
www.punktal.com